


Evolution. Revolution.

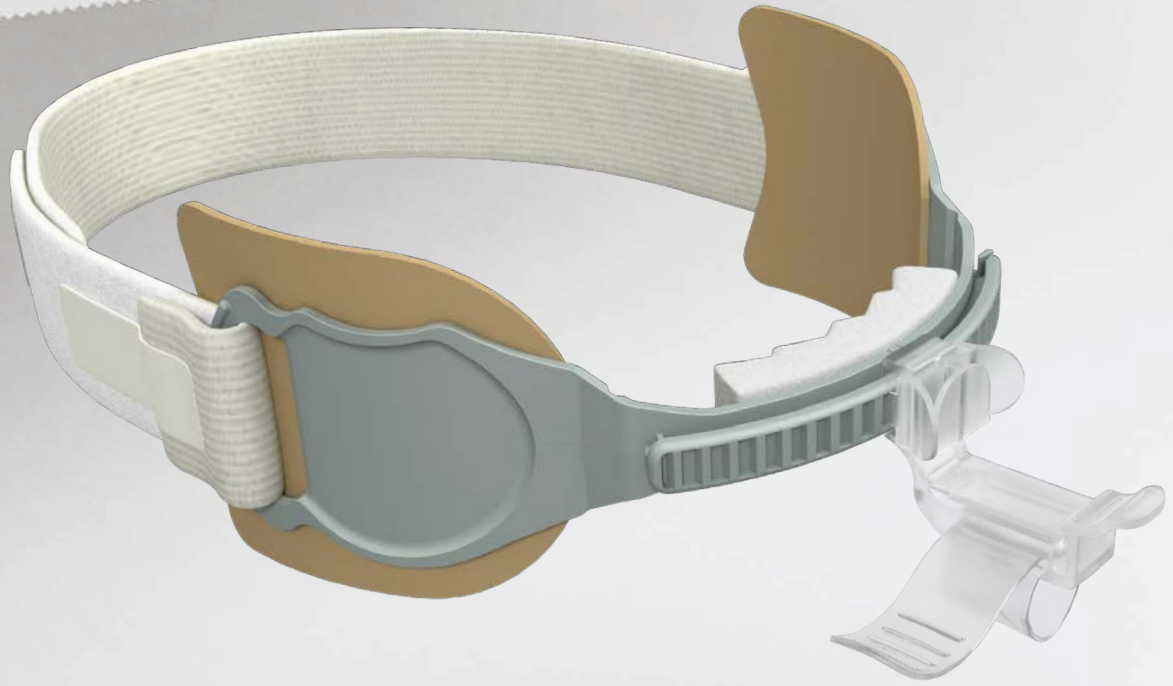

Innovated for you — to leave tape behind.


The AnchorFast SlimFit™ oral endotracheal tube fastener provides a convenient means to hold an oral endotracheal (ET) tube securely in place without the use of adhesive tape and is designed with the intent to accommodate smaller faces.

AnchorFast SlimFit
Oral Endotracheal Tube Fastener

AnchorFast SlimFit

Oral Endotracheal Tube Fastener

Designed to accommodate smaller faces

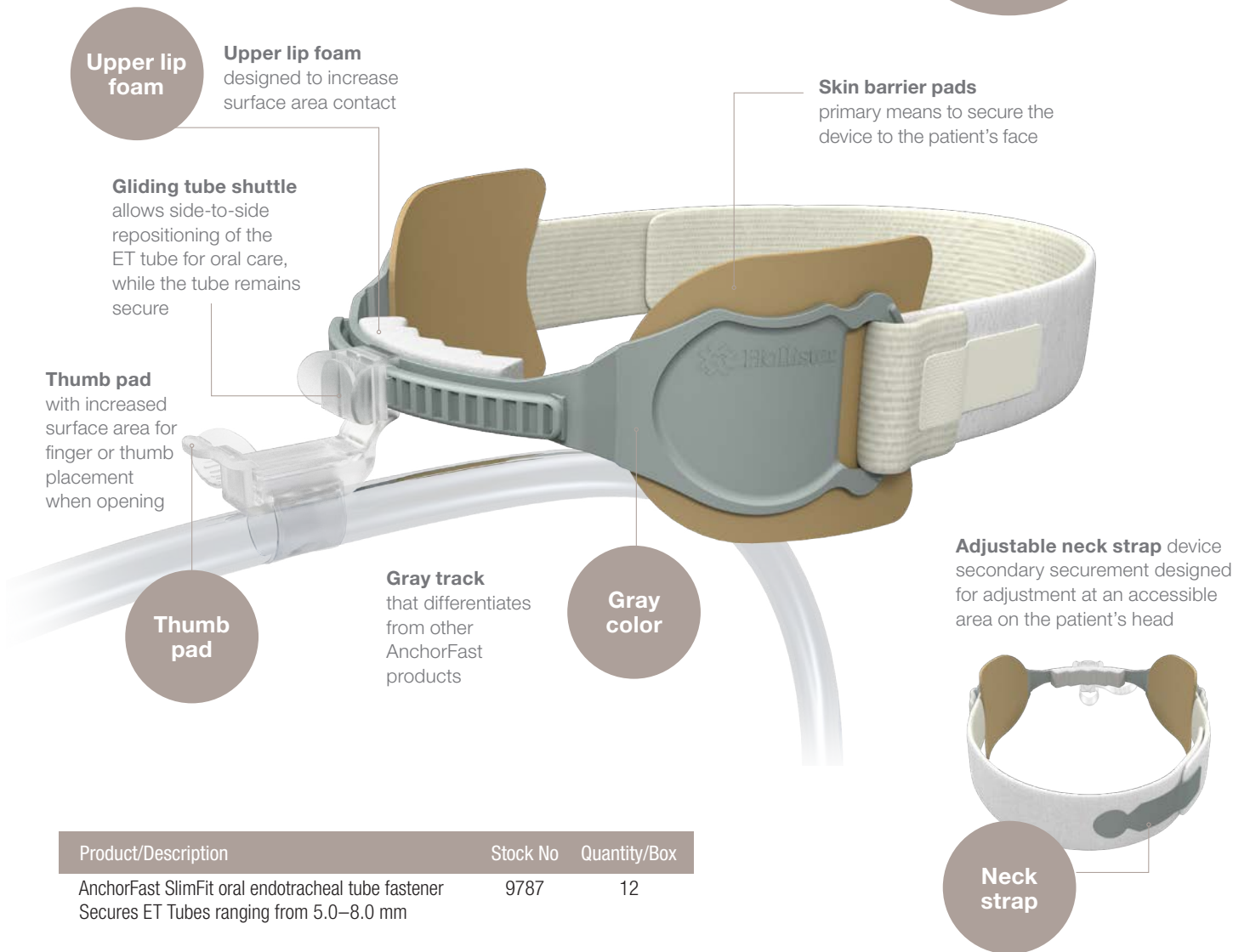

Upper lip foam

Upper lip foam designed to increase surface area contact

Gliding tube shuttle allows side-to-side repositioning of the ET tube for oral care, while the tube remains secure

Thumb pad with increased surface area for finger or thumb placement when opening

Thumb pad

Gray track that differentiates from other AnchorFast products

Gray color

Skin barrier pads primary means to secure the device to the patient's face

Adjustable neck strap device secondary securement designed for adjustment at an accessible area on the patient's head

Neck strap

Product/Description	Stock No	Quantity/Box
AnchorFast SlimFit oral endotracheal tube fastener Secures ET Tubes ranging from 5.0–8.0 mm	9787	12

For more information or to request a sample, visit hollister.com/anchorfast



Hollister Incorporated
2000 Hollister Drive
Libertyville, Illinois 60048
1.800.323.4060

www.hollister.com

USA Rx Only. Caution: Federal (USA) Law restricts this device to sale by or on the order of a physician or other healthcare practitioner licensed under state law to order this product. Prior to use, be sure to read the Instructions for Use for information regarding Intended Use, Contraindications, Warnings, Precautions, and Instructions.

The Hollister logo and AnchorFast SlimFit are trademarks of Hollister Incorporated. © 2020 Hollister Incorporated.

923320-1020

